


# CompreDrive Distance add-on for CompreCell Pouch User Manual

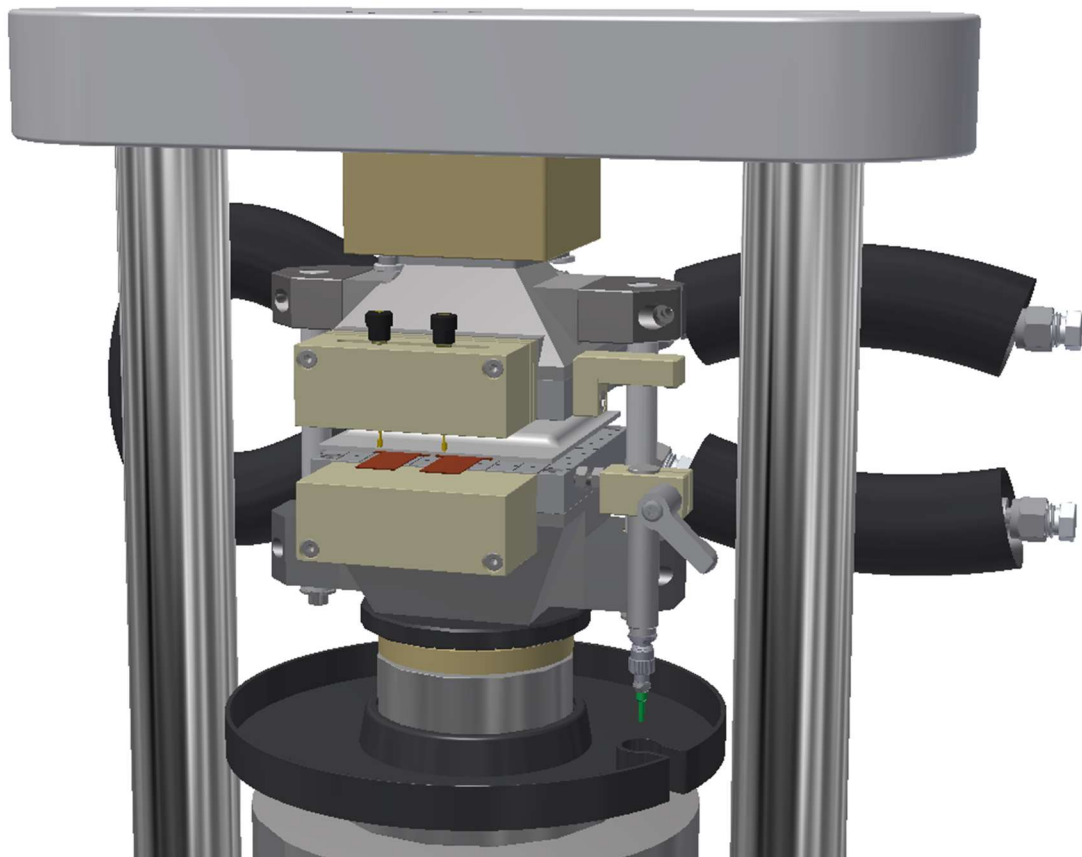

Version 1.2 EN

03/2024

rhd instruments GmbH & Co. KG



# Contents

1	Product description	1
2	General information	1
3	Important general safety notes	3
4	Components of the CompreDrive Distance add-on system	5
5	Operation conditions, storage and rated values	7
6	Essential features at a glance	8
7	Getting started	9
	7.1 Retrofitting the CompreDrive Distance add on	9
	7.2 Replacement of the lower drip pan with the Distance Add-On version	13
	7.3 Sensor and mirror mount assembly	13
8	Rough adjustment of the sensor height	19
9	Initial software setup	22
10	Fine adjustment of the sensor height	27
11	Troubleshooting	31
	11.1 No Connection to confocal system, system not found	31
	11.2 Confocal system connected, but no data is transferred to CompreDriveControl	31
	11.3 Data is present at the start point of a measurement, but is missing from some point	31
12	Settlement	32
13	Contact and Technical Support	33

# 1 Product description

The CompreDrive Distance add-on is designed to enable the electrochemical characterization of pouch cells while precisely monitoring the thickness changes of the sample during measurement. The product is designed to be used only in combination with the corresponding measuring cells (CompreCell Pouch) and the CompreDrive setup.

The CompreDrive Distance add-on provides an optical confocal distance sensor which can be directly mounted on the CompreCell Pouch. The small and light-weight optical sensor is connected via an optical cable to the confocal analysis unit. This unit can be directly attached to the CompreDrive main frame, where power is supplied. The data is transferred via an ethernet interface and can directly be recorded alongside with all other data using CompreDriveControl software.

The optical sensor offers a total measuring range of 2 mm with a resolution of 40 nm and a linearity of 1  $\mu\text{m}$  over the whole range.

## 2 General information



Thank you for your confidence in our products and services. We wish you pleasure and success with your new CompreDrive Distance add-on system. It has been especially been developed for professional electrochemical material characterization.

- » **To avoid physical injuries and damages, please read this instruction manual carefully before using the device for the first time.**
- » **Please pay attention to all safety notes in this instruction manual.**
- » **Please keep this manual safe. In case of selling or leaving the device to third parties, please do not forget to hand this manual over as well.**
- » **The operation of the CompreDrive Distance add-on should only be performed by properly trained and experienced members of staff.**
- » **The setup is developed to measure sample thickness changes in your CompreCell electrochemical measuring cells and must not be used for any other purpose.**

- » To avoid unstable operating conditions and injury, the CompreDrive Distance setup as well as the individual components should not be used if
- they show noticeable damage,
  - they were stored or operated under unapproved conditions (see operational condition, storage and rated values),
  - they were exposed to high mechanical stress, exceeding normal usage,
  - they were altered by members of staff not authorized by rhd instruments.

The instructions in this manual were carefully checked for correctness. However, liability for any mistakes in form and content will not be assumed. Additionally, rhd instruments GmbH & Co. KG (in the following declared as rhd instruments) reserves the right to change the setup and design of the products presented and described within this manual. Such changes are necessary to guarantee the continuous development of the products and, thus, the improvement of product quality and reliability.

### Markings in this manual


Marking	Meaning
 <b>WARNING</b>	Indicates a hazardous situation which, if not avoided, could result in a serious injury or death.
 <b>ADVICE</b>	Indicates potential physical damages and other important information associated with your device.

### 3 Important general safety notes



**WARNING: Do not bend or kink optical fibres**

If the optical cable between sensor and controller is kinked or bended in tight radii, there will be irreversible damage to the cable. Minimum radius is 40 mm.

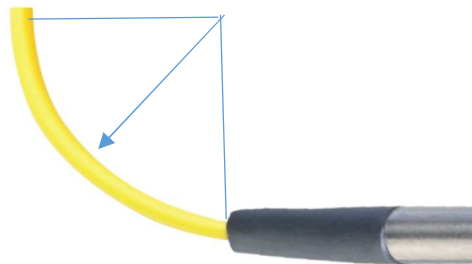

- » Connect the power supply and the display/output device according to the safety regulations for electrical equipment. Otherwise, there is risk of injury, damage to or destruction of the sensor and/or the controller
- » The supply voltage must not exceed the specified limits to avoid damage to or destruction of the sensor and/or the controller
- » Protect the ends of the optical fibre against contamination (use protective caps).
- » Only use the original parts included in delivery. They are prepared for your device and guarantee the necessary safety for operator and device.
- » Do not operate the device with wet hands. Operate the device only in dry rooms.
- » Do not operate the device outdoors.
- » Please follow only the instructions in this manual for cleaning the device.
- » Make sure that cables and conductors are not damaged. Damage could be caused by heat, impact, contact with chemicals, or mechanical impacts like rubbing, bending, tearing, and rolling-over.

- » Prevent the device from mechanical impact. In case the device fell down, please contact rhd instruments or a technician authorized by rhd instruments before switching it on again.
- » If your device shows any visible damage or defect: Disconnect the power supply by pulling out the power connector. Never operate your device in a damaged state. Never repair the device on your own. The device should only be repaired by either rhd instruments or by a technician authorized by rhd instruments.
- » Do not open the device. There are no user-serviceable parts inside.
- » Please follow this instruction manual for maintaining your device.
- » Only use original spare parts delivered by rhd instruments.



**ADVICE: Connect and disconnect any cable connection carefully.**



**ADVICE: Handle chemicals with care.**

- » When handling chemicals during preparation and execution of measurements with the CompreDrive HC add-on, the usual safety advice in accordance with the H, EUH, and P statements (in the European Union: rating principles according to the CLP regulation) and appropriate safety measures have to be observed. This applies to subsequent cleaning and decontamination as well.



**ADVICE: Sufficient cleaning increases the lifetime of your system.**

- » After using the CompreDrive add-on for electrochemical measurements, all components in contact with chemicals need to be thoroughly cleaned. Insufficient cleaning, decontamination, and drying of the components may result in damage due to corrosion and, thus, may affect the quality of your measurement results.

## 4 Components of the CompreDrive Distance add-on system

- » Please unpack your device carefully.
- » Please check if the delivery is complete:
  - 1x confocal controller unit with support plate
  - 1x optical sensor with optical cable and holder
  - 1x power cable 24 V (Y-cable)
  - 1x ethernet cable 2 m
  - 1x hex screw driver (4 mm)
  - 1x wrench 8 mm
  - 1x sensor mount
  - 1x mirror mount
  - 1x wedge
  - 1x manual "CompreDrive Distance add-on"
- » Please check if the delivered items are undamaged.



**ADVICE:** If the delivered items are incomplete or damaged, please contact rhd instruments via e-mail ([info@rhd-instruments.de](mailto:info@rhd-instruments.de)) or via phone (+49-6151-8707187).

rhd instruments will reject any claims for warranty or responsibility in case damaged equipment is used.

In case any accessory of other manufacturers is used, rhd instruments will accept no liability.





Figure 1: Overview CompreDrive Distance add-on for CompreCell Pouch.

# 5 Operation conditions, storage and rated values

- » Power supply 24 VDC  
 $\pm 15\%$ ,  $I_{max} < 1\text{ A}$
  
- » Temperature range during operation  
 Sensor:  $T_{env. operation} = -20\text{ }^{\circ}\text{C}$  to  $+70\text{ }^{\circ}\text{C}$
  
- » Temperature range during operation  
 Controller:  $T_{env. operation} = +5\text{ }^{\circ}\text{C}$  to  $+50\text{ }^{\circ}\text{C}$
  
- » Temperature range during storage:  $T_{storage} = +10\text{ }^{\circ}\text{C}$  to  $+50\text{ }^{\circ}\text{C}$
  
- » Relative humidity (RH) for working (non-condensing) 0 to 80% RH  
 and storage area:
- » Atmosphere during storage: Non-corrosive


**Confocal sensor characteristics**

Measuring range	2 mm
Start of measuring range approx.	14 mm
Lateral resolution*	40 nm
Light spot diameter	10 $\mu\text{m}$
Linearity (displacement and distance measurement)	1 $\mu\text{m}$

\* Average from 512 values at 1kHz, near to the midrange

## 6 Essential features at a glance

- » High-precision distance measuring system for CompreCell Pouch.
- » Measurement of sample thickness changes during electrochemical experiments.
- » Measuring points as close to the sample as technically possible, to exclude influences from the main frame.
- » Fast and comfortable assembly and handling.
- » Fully included data handling and processing into CompreDriveControl software.
- » Retrofittable to all CompreDrive systems.
- » Compatible with HC temperature control option.

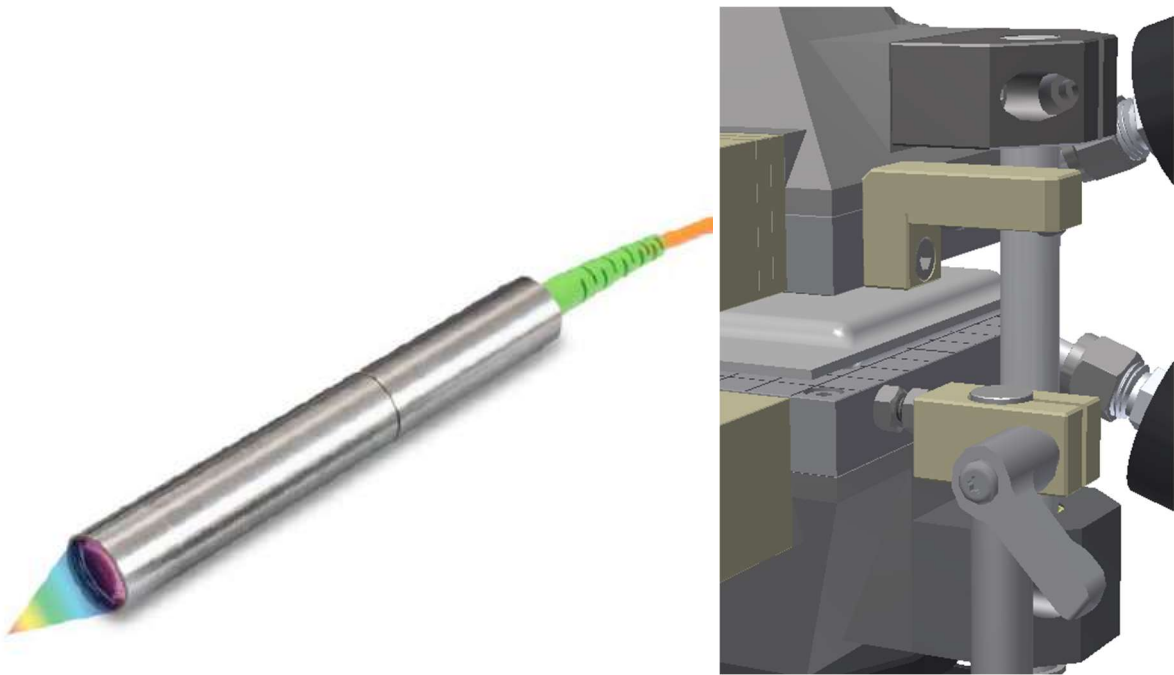

### Note:

If you have any questions, for example with regard to the compatibility of your measurement devices, do not hesitate to contact us via email ([info@rhd-instruments.de](mailto:info@rhd-instruments.de)) or phone (+49-6151-8707187).

## 7 Getting started

### 7.1 Retrofitting the CompreDrive Distance add on


If your CompreDrive was not delivered with the CompreDrive Distance option by default, it needs to be installed as described in the following:

**Switch the CompreDrive off at the main switch!**



**Figure 2: Confocal controller on support plate with aligned sliding blocks.**

If the confocal controller is not yet mounted on the support plate, click it onto the top hat rail. Then make sure, that the pre-installed sliding blocks are aligned as shown.

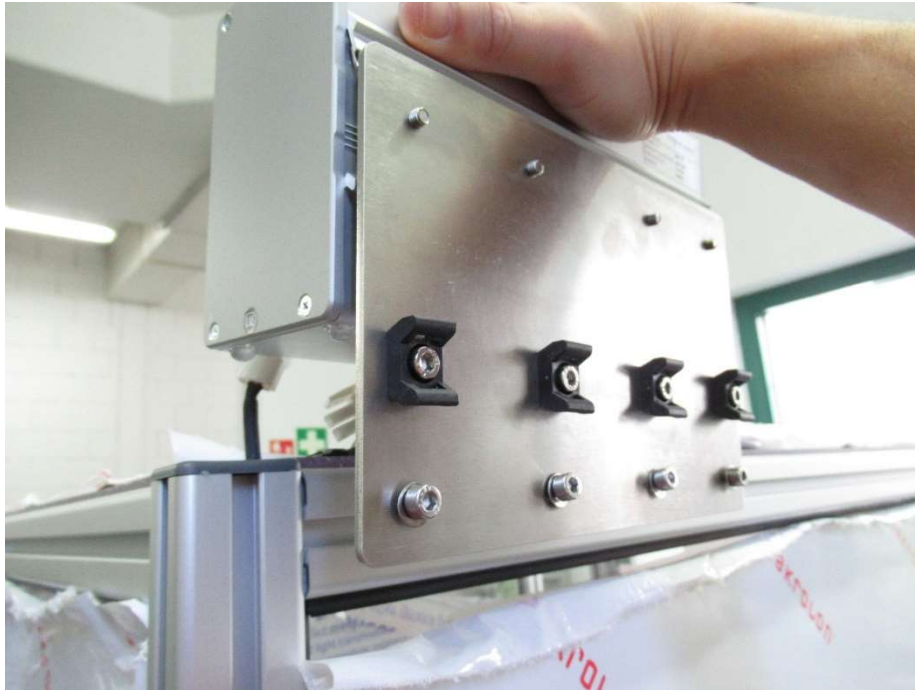

**Figure 3: Mounting position of confocal controller**

The unit is fixed at the top left corner of the back-side of the CompreDrive. Push the sliding blocks into the aluminium profile and use a 4 mm hex screw driver to fix the plate.



**Figure 4: Confocal controller mounted on the CompreDrive**



**Figure 5: Left: Y-cable connected to 24 supply and to the fan, right: Connection to controller**



**Figure 6: Sensor cable inserted into the sensor port of the controller**


Put the optical cable through the opening in the back panel. Then remove the black protecting plug from the sensor port of the controller and insert the sensor connector. It must fully engage.



**Figure 7: Fixing points for all cables at the backside of the sensor support plate.**

Now, all cables can be fixed to the back side of the sensor support plate. Use cable ties.


## 7.2 Replacement of the lower drip pan with the Distance Add-On version

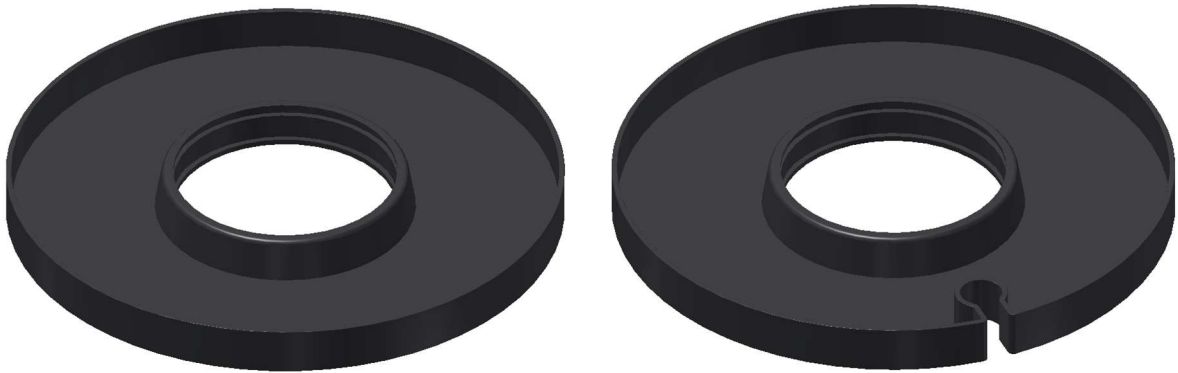

**Figure 8: Standard drip pan vs. drip pan for Distance Add-On**

For instructions on how to remove and place the drip pan see CompreCell Pouch 10S HC manual.


## 7.3 Sensor and mirror mount assembly

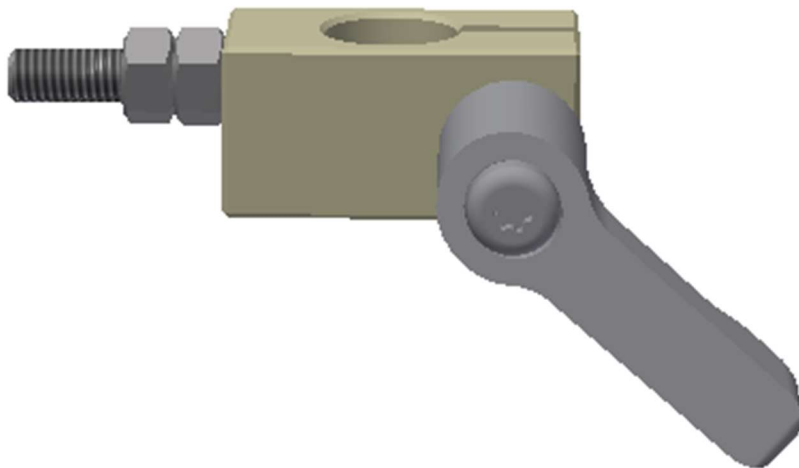

**Figure 9: Sensor mount**


Make sure that the two nuts are screwed towards the PEEK part of the sensor mount as far as possible.



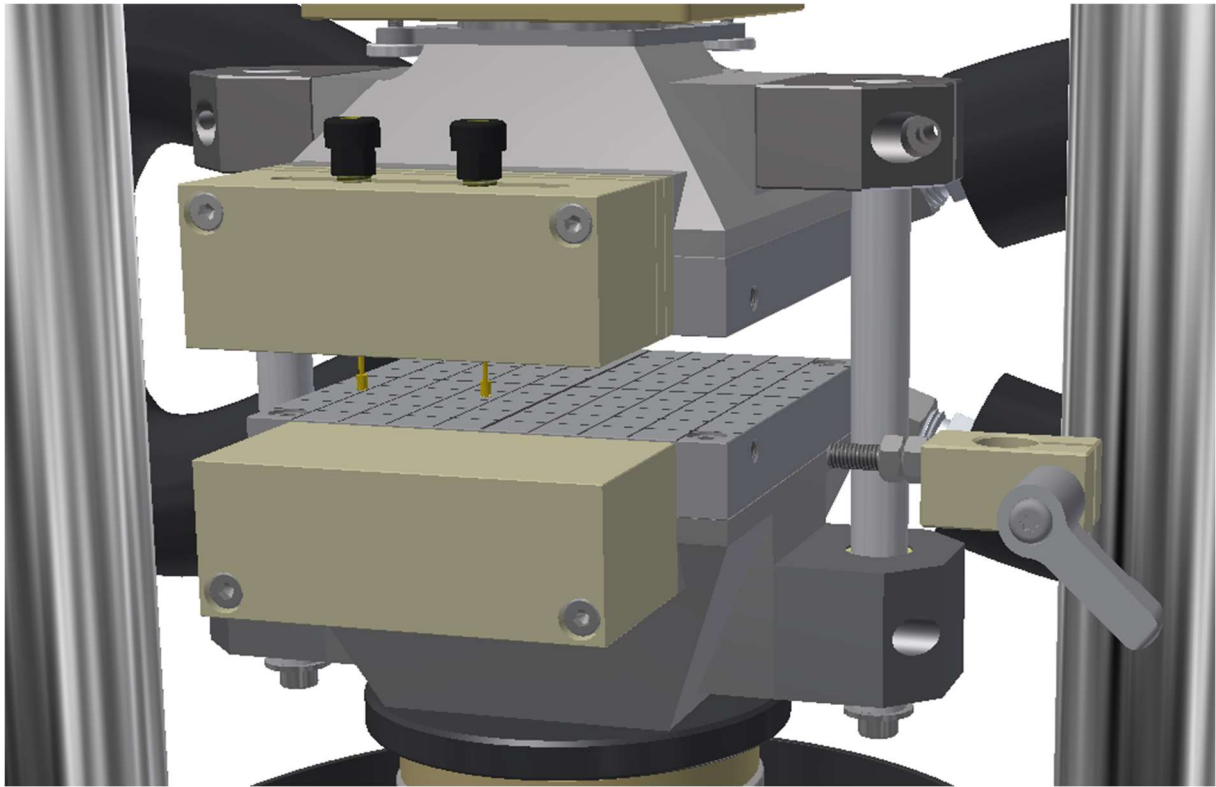


**Figure 10: Attaching the sensor mount to the CompreCell Pouch**

Attach the sensor mount to the CompreCell Pouch by screwing it into the lower, right screw hole shown in **Figure 10**.

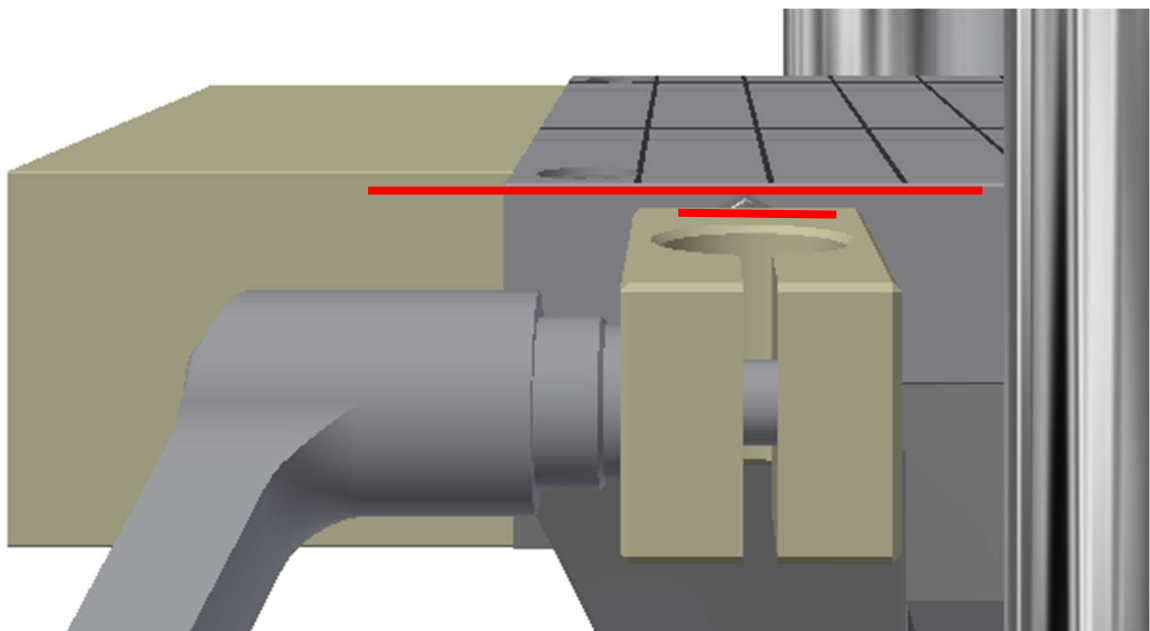


**Figure 11: Aligning the sensor mount**

Screw the sensor mount in as far as possible, then turn it back so that the upper surface of the PEEK part is parallel to the ground surface on your CompreCell Pouch.

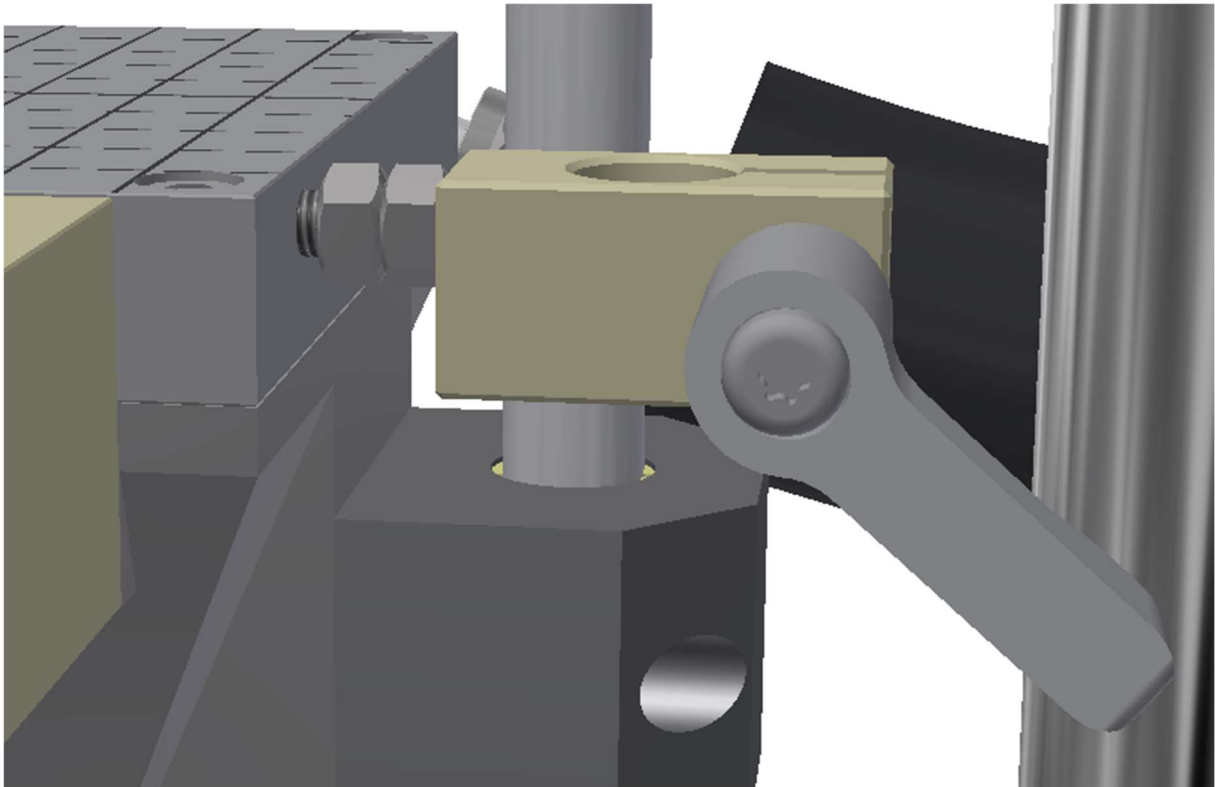


**Figure 12: Nuts to be tightened**

Lock the sensor mount in place by screwing the nut closer to your CompreCell Pouch against the corresponding surface with an 8 mm wrench while holding the sensor mount in its desired position.

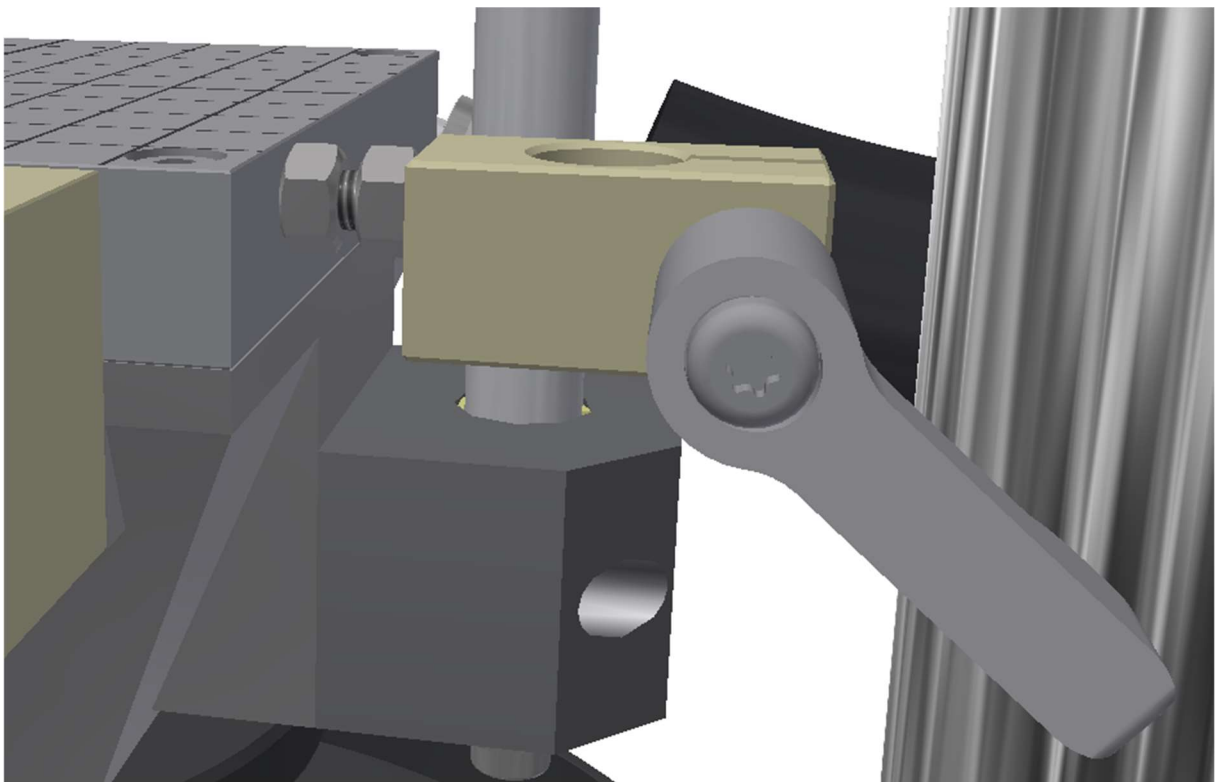


**Figure 13: locked sensor mount**

The sensor mount should now be locked in place securely.

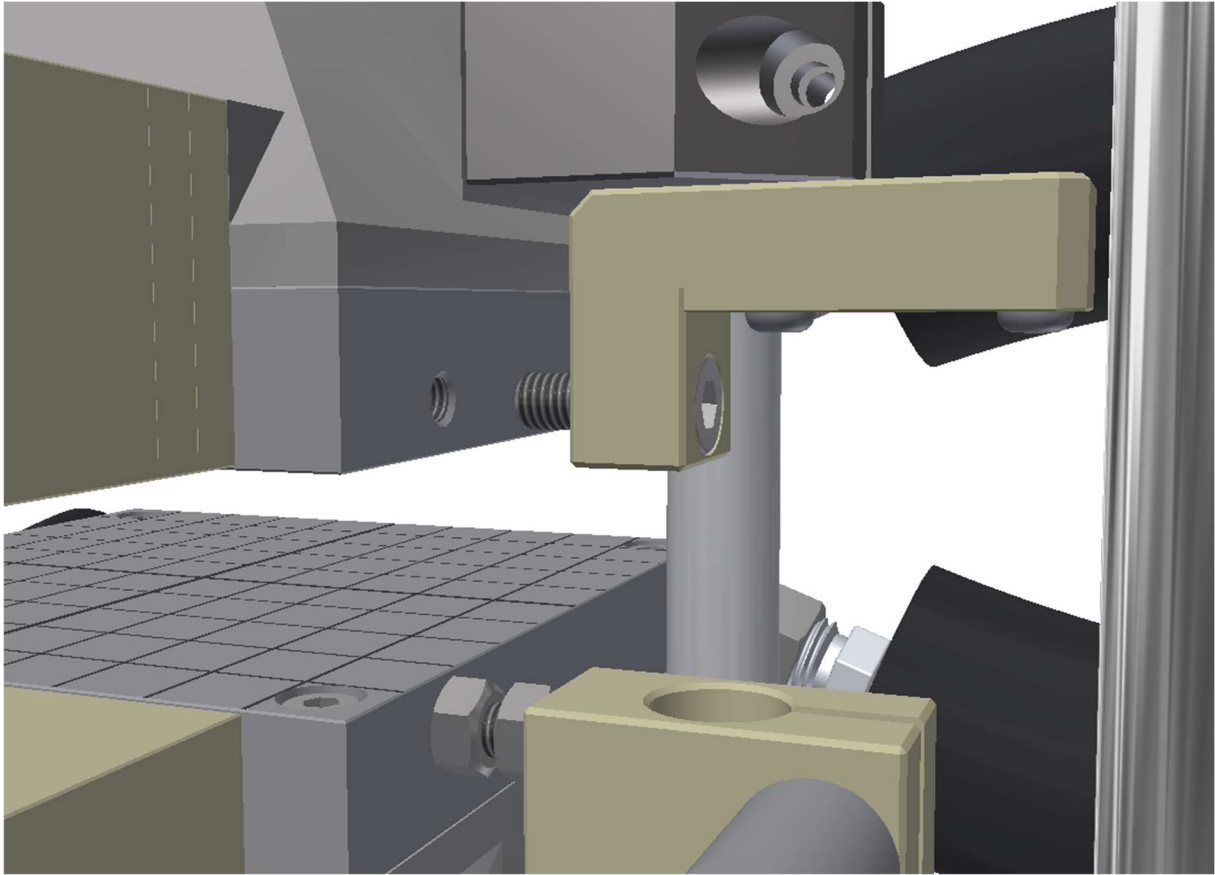


**Figure 14: Attaching the mirror mount**

Attach the mirror mount to your CompreCell Pouch with the included M5 screw and a 4 mm hex key.

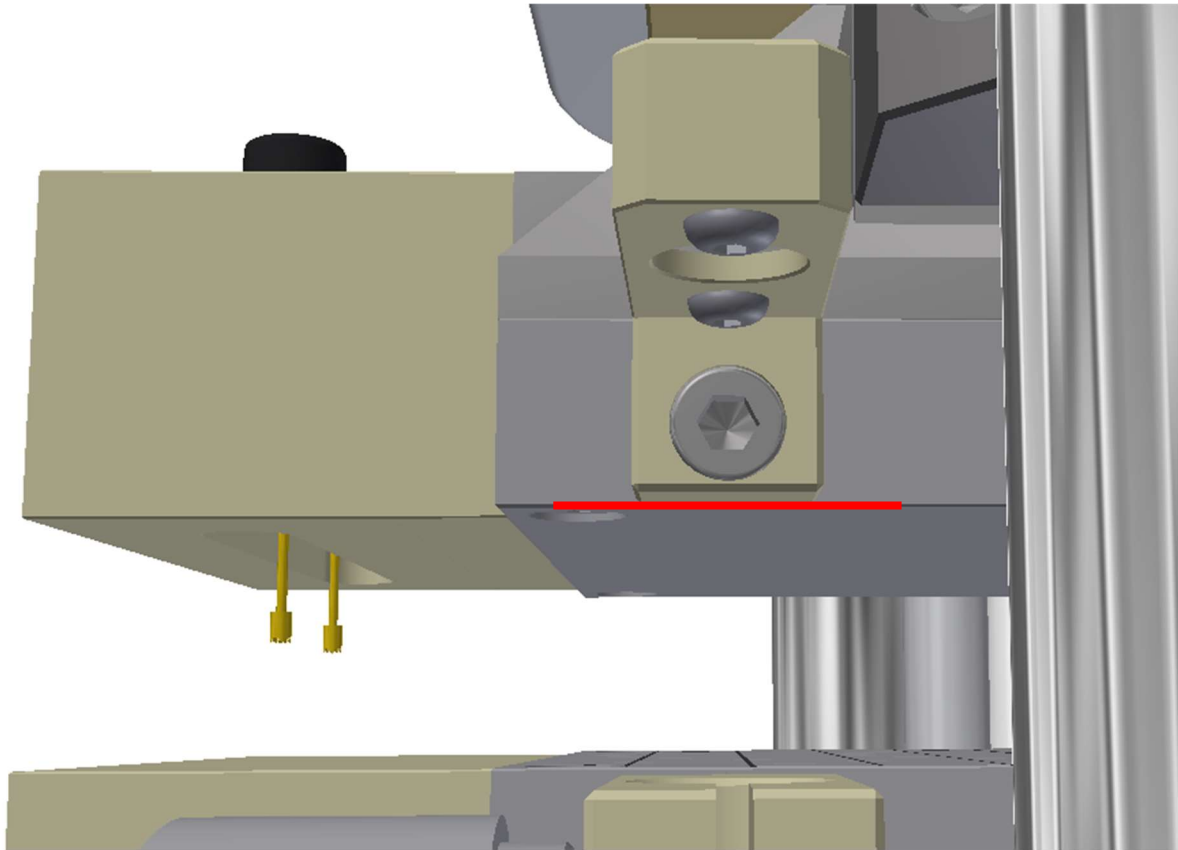


**Figure 15: Aligning the sensor mount**

Make sure that the bottom part of the sensor mount is parallel to the upper plate of the CompreCell Pouch. It must be even with the surface of the upper plate, not protruding.

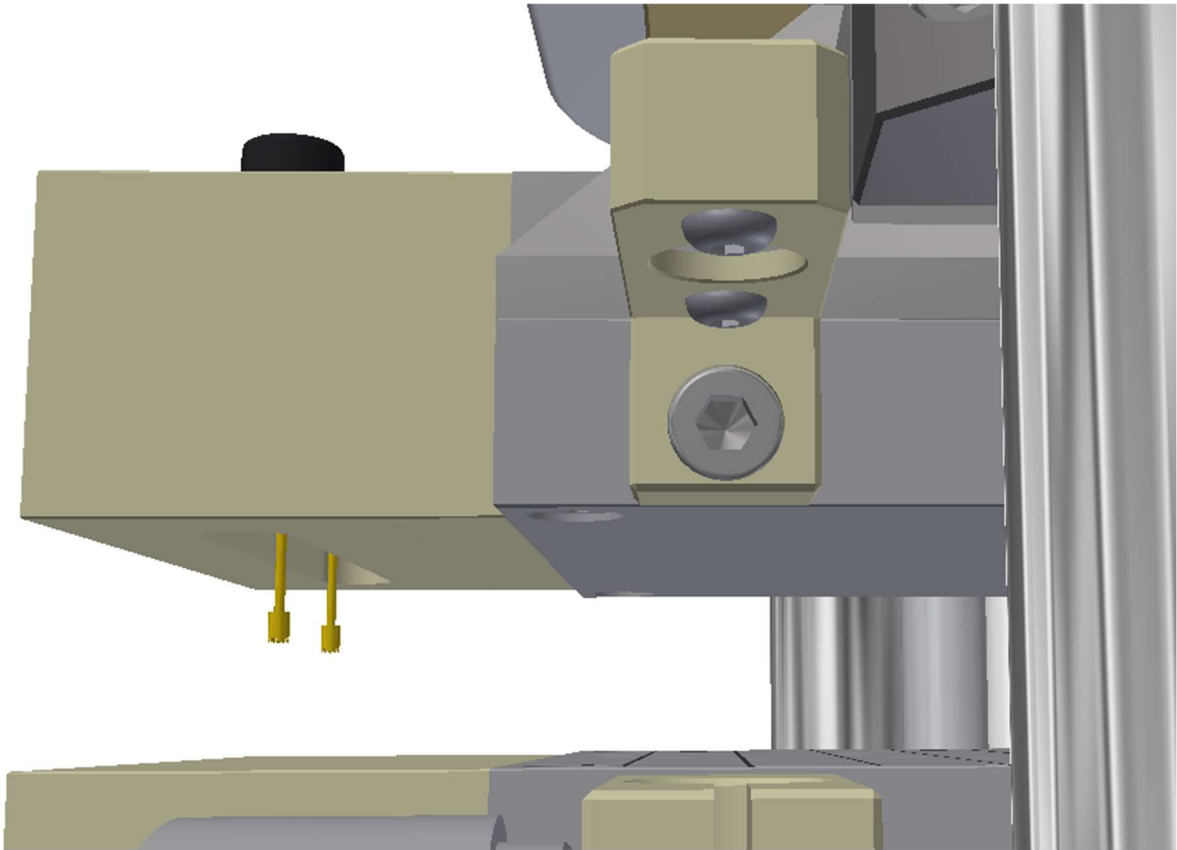


**Figure 15: Fixed sensor and mirror mounts**

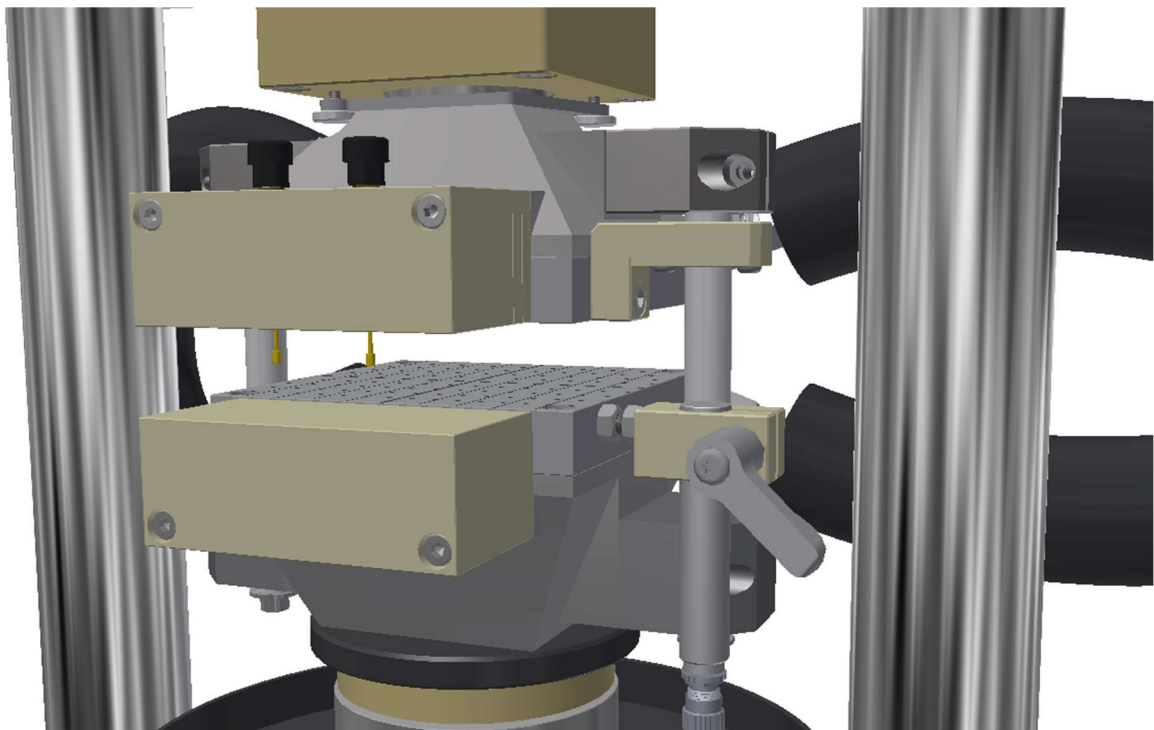

**Figure 16: Inserting the confocal sensor**


Insert the confocal sensor into the sensor mount and lock it in place by turning the lever on the sensor mount.

## 8 Rough adjustment of the sensor height

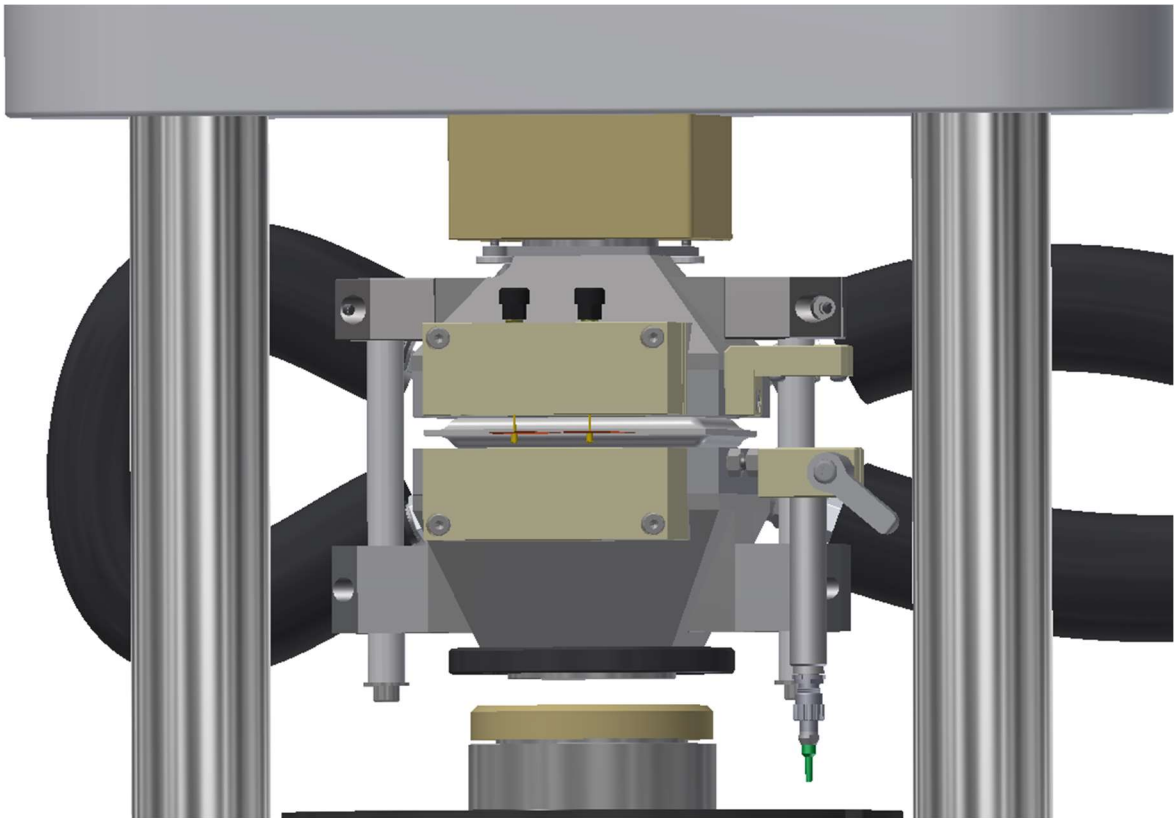

**Figure 17: Distance add-on mounted completely**

The distance between sensor and mirror has to be adjusted according to the height of the pouch cell in use. To do so, start by lifting the bottom part of your CompreCell Pouch until pouch battery cell is in full contact.

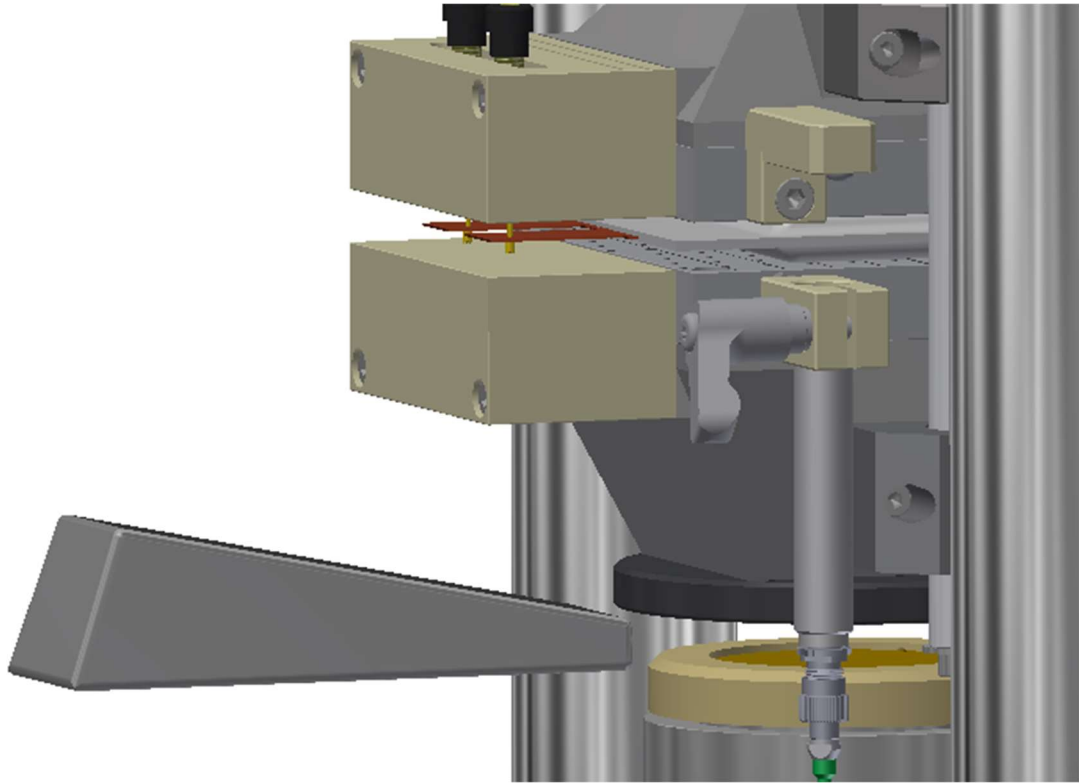

**Figure 18: Inserting a wedge for easy sensor adjustment**

Place the wedge between the PEEK ring and the bottom part of the CompreCell Pouch to hold the CompreCell Pouch in position.



**Figure 19: Sensor spot visible on the mirror.**

Move the sensor until a small sharp light spot becomes visible on the mirror.

**Hint:** The focal distance of the confocal sensor is 14 mm. Since the mirror is recessed about 2 mm in the mirror mount, the distance between the sensor and the mirror mount should be roughly 12 mm.

Make sure that both sensor LEDs on the confocal sensor controller are green at the same time.



**Figure 20: Sensor LEDs on the sensor controller.**

Fine tuning of the sensor position and measuring range is done using the RHD.MECServer software, see next chapters.

**As described in chapter 10, make sure to remove the wedge from below your Compre-Cell Pouch before approaching the CompreDrive!**



## 9 Initial software setup

The CompreDriveControl software can be used to establish automated electrochemical measurements including the communication with the confocal sensor system. For further details please refer to the CompreDriveControl user manual as well as the manual of the confocal controller.

To establish a connection, the RHD.MECServer plugin has to be installed. Unpack the RHD.MECserver archive. Copy the files from the "Plugin" folder of the archive to the following direction:

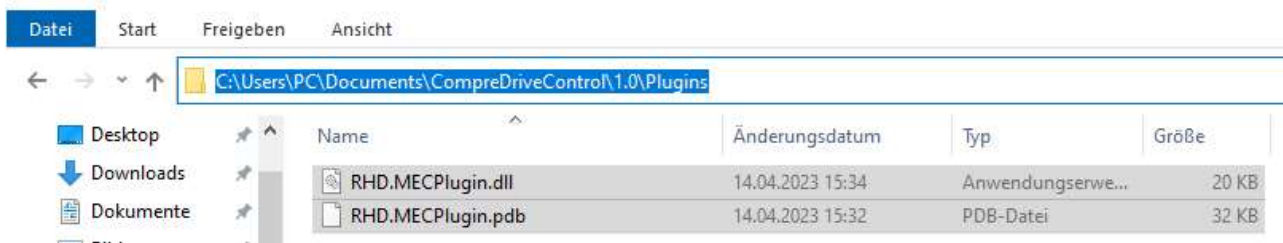


Figure 21: Files to be copied to the Plugin folder of CompreDriveControl.

New plugins are only recognized at start-up of CompreDriveControl.


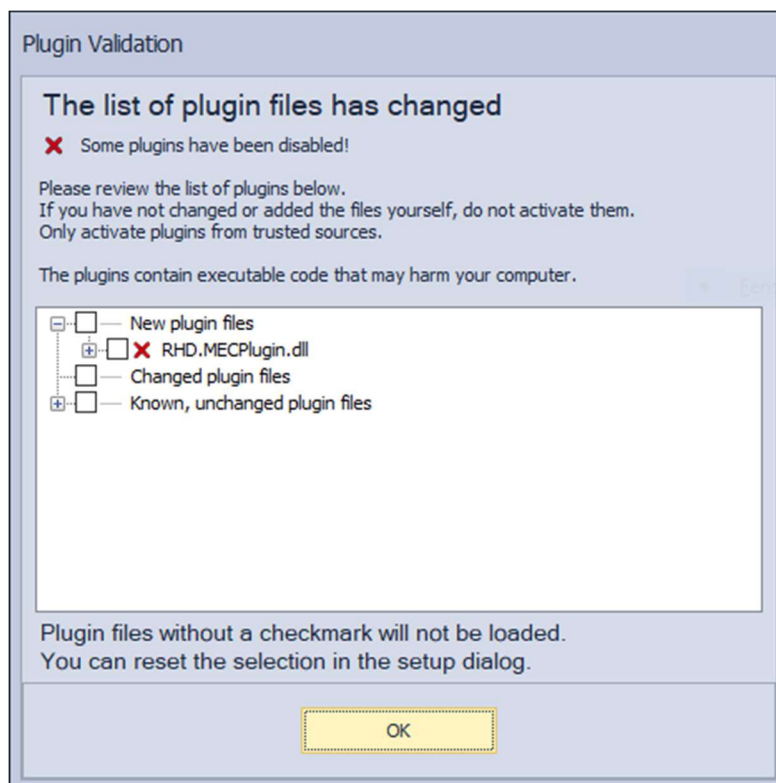



Figure 22: New plugin recognized. Click the checkbox and allow it to be enabled

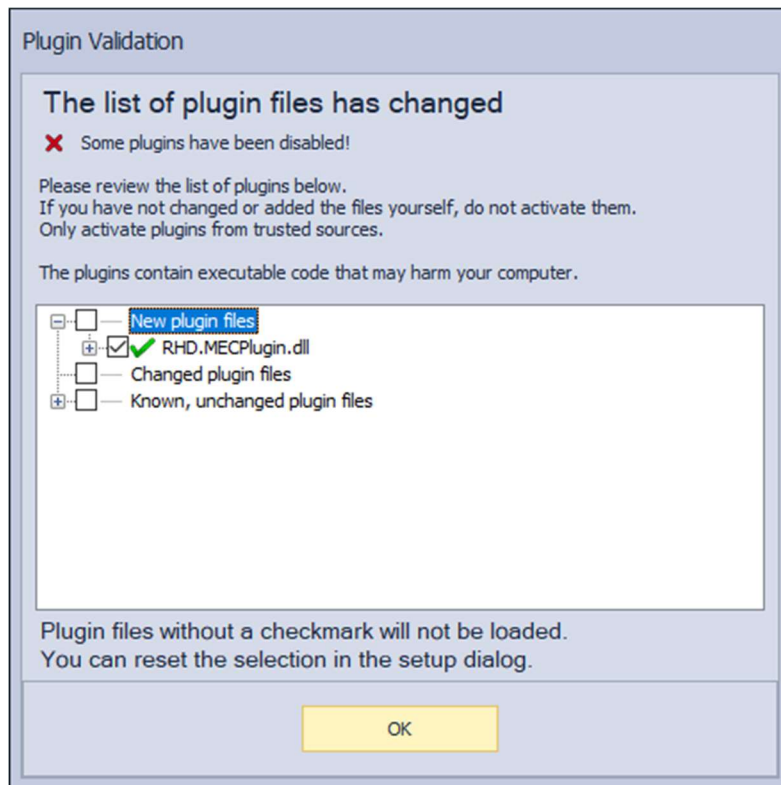

**Figure 23: New plugin enabled**


To make sure the plugin is recognized, open CompreDriveControl go to "Setup" in the Main Ribbon, and then click on "Control Plugins".

In the list should all available data provider plugins be visible. "RHD MEC addon RAW" and "RHD MEC addon TARA" must be visible and activated.

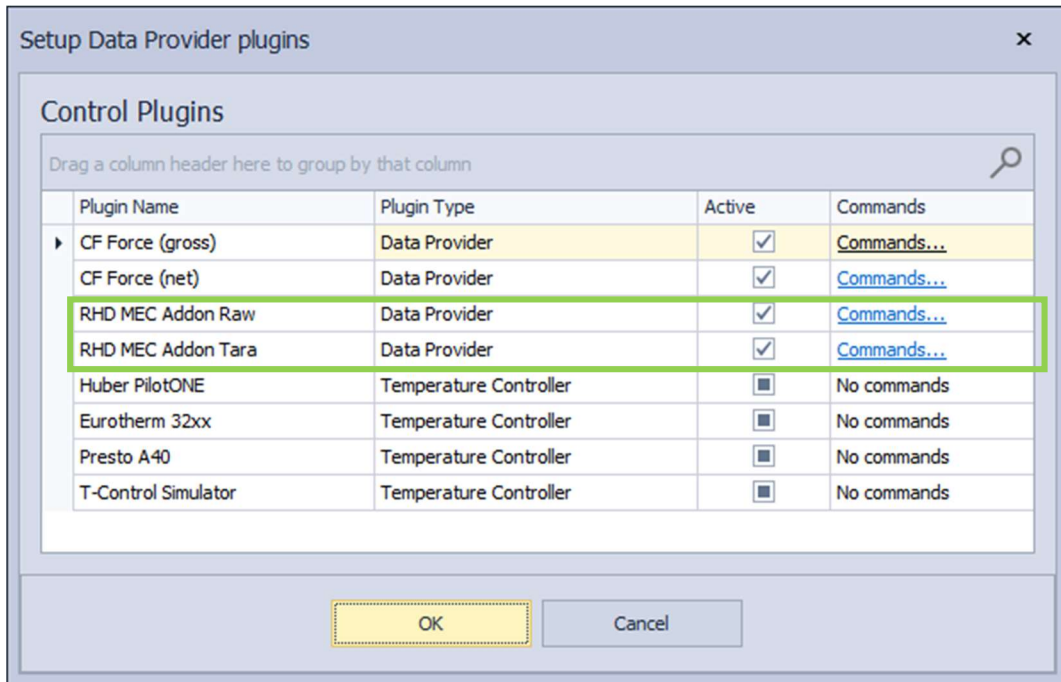


**Figure 24: Selection of data providers, confocal distance is recognized.**

Close CompreDriveControl and start the MECServer software by executing RHD.MECServer.exe.

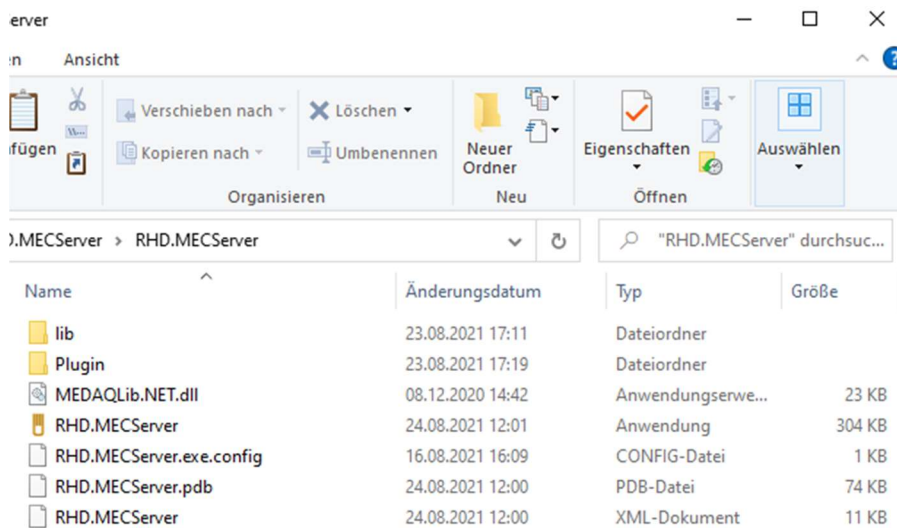

**Figure 25: Content of the RHD.MECServer archive.**

The following window opens:

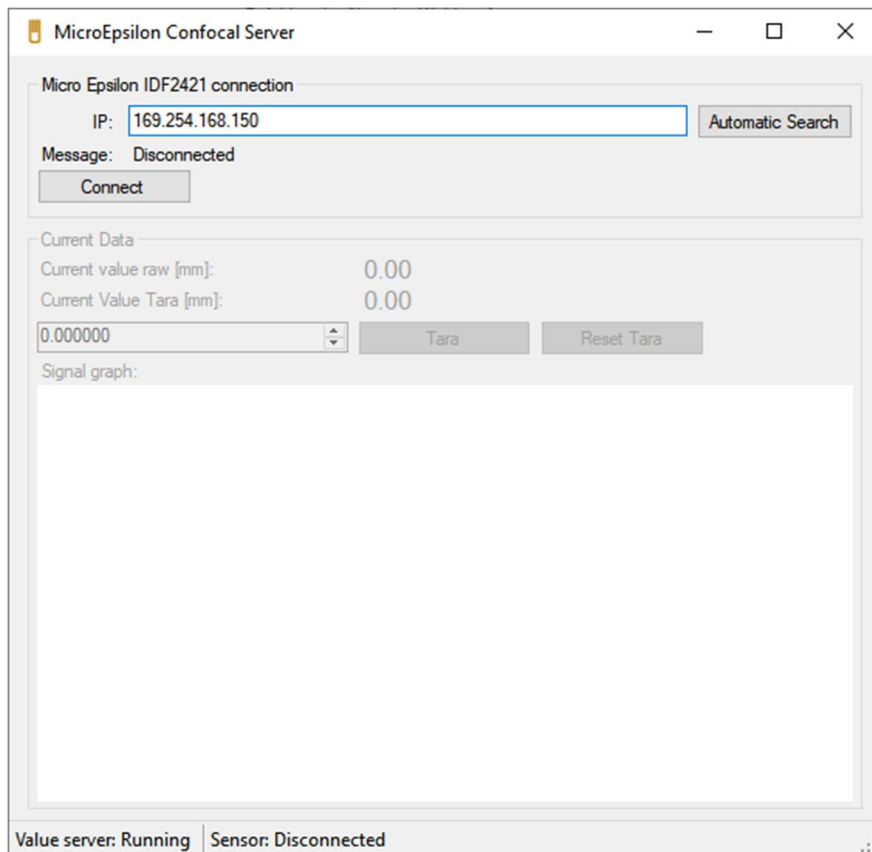


**Figure 26: Startup window of RHD.MECServer.exe**

The default IP address of the controller is **169.254.168.150** (static). It is recommended to use the provided USB-to-ethernet adapter to create a direct connection between the controller and the PC. The adapter should already be configured to use the same subnet as the controller.

In case a connection to the default IP address does not work, click the **Automatic Search** button. Here, the software will automatically scan the current subnet for controllers. If the controller can't be found, please check the following steps:

- Connect the sensor directly via the supplied USB-Ethernet dongle
- Ensure in the network settings, that the USB-Ethernet dongle has a static IP of 169.254.168.10 assigned
- Ensure that TCP/IP communication via port 1024 is not blocked by a firewall
- After plugging in the cable, make sure the two LEDs on the controller connector light up (green and yellow)  
→ If not, it means that the connection is not established correctly. Check the cable for damage.

After entering the IP address of the controller, click on "Connect". If the controller is switched on, connected to the network, and the right IP is set, connection is established. As soon as the connection works, data from the controller is transferred and displayed.

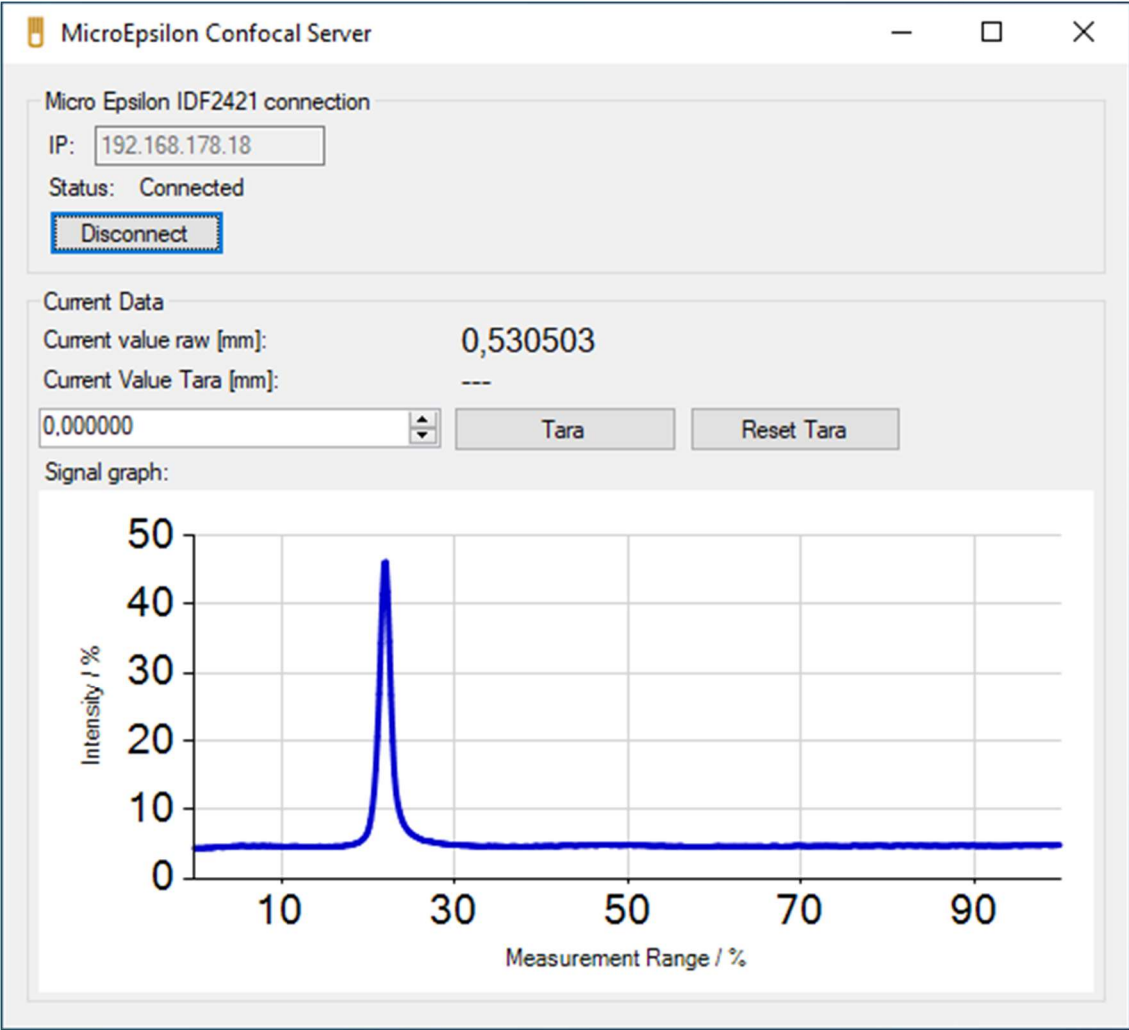

**Figure 27: Confocal Server connected to confocal controller.**

Start CompreDriveControl, data from the confocal system should now be available.

# 10 Fine adjustment of the sensor height


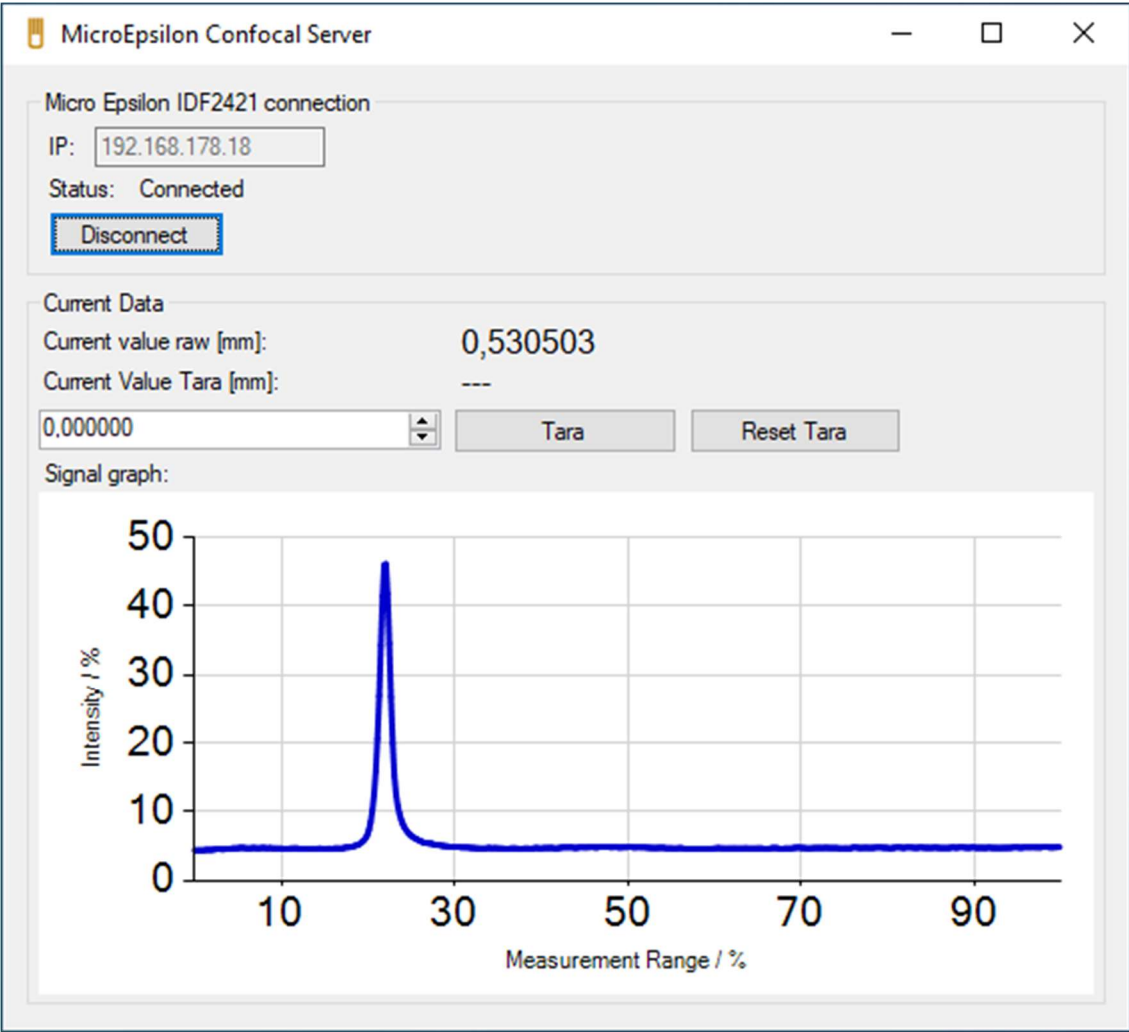



Figure 28: Desired sensor signal

A single sharp peak is visible in a range between 10 % and 90 % of the measurement range is visible, if the sensor is within its working limits (Figure 28). This is also indicated by two green LEDs on the controller panel.

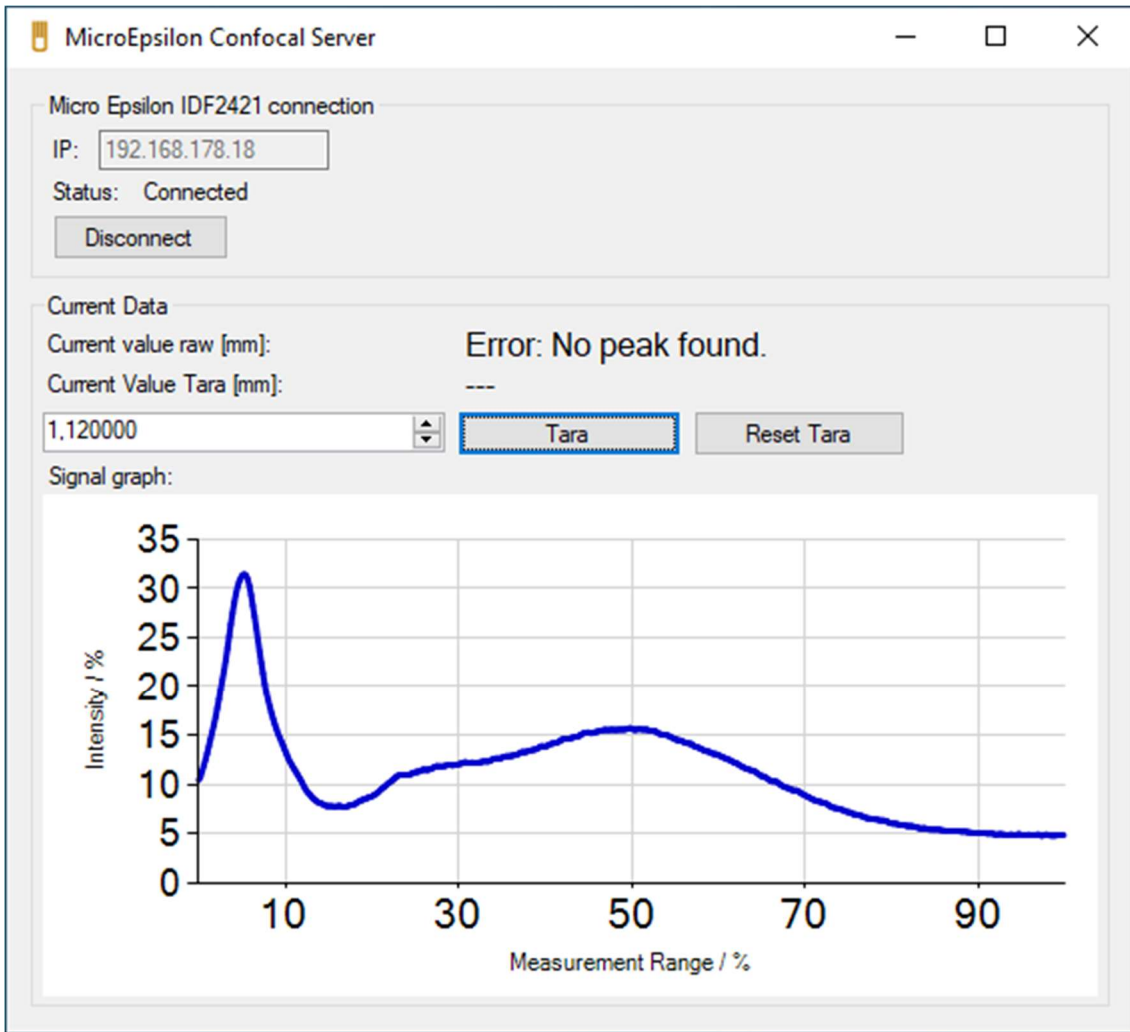

**Figure 29: Sensor signal off-limit**

If the sensor is off-limit (protective cap in place, out of distance), a result like in the above picture is visible (**Figure 29**).

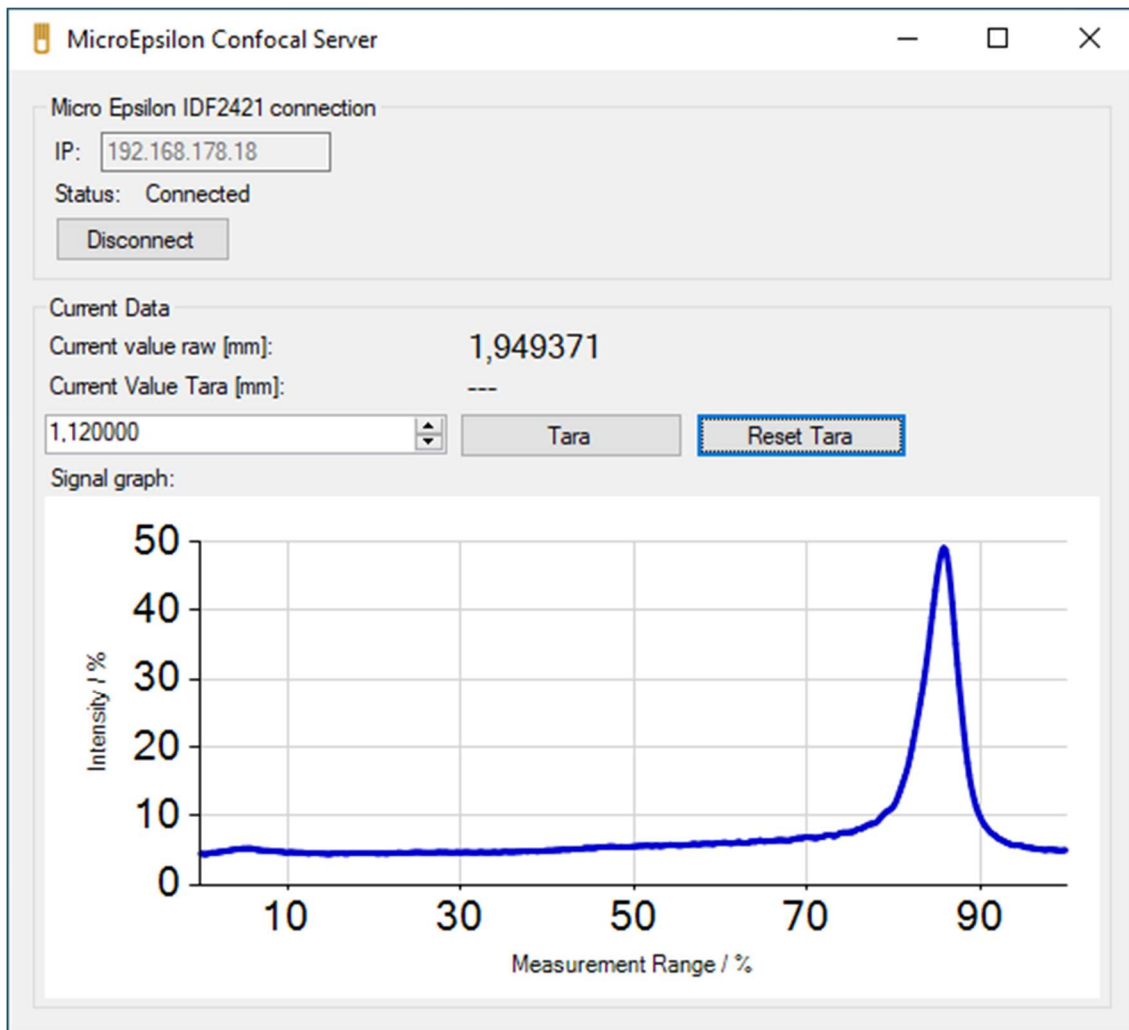

**Figure 30: Setting the sensor to the uppermost limit to use full range of compression.**

If a highly compressive sample is loaded, make sure to use the whole range by adjusting the sensor to the uppermost limit prior to measurement, as shown in the picture below (**Figure 30**).

If the initial thickness of the sample is known from comparative measurements of the empty cell vs. the loaded cell (using a thickness measuring gauge on a plane plate), tara can also be set.

Lift the bottom part of your CompreCell Pouch to lift the weight off of the inserted wedge. **REMOVE THE WEDGE BENEATH YOUR COMPRECELL POUCH!**

**WARNING: Leaving the wedge in between the CompreDrive and the CompreCell Pouch while using the CompreDrive can cause permanent damage to the equipment.**

When using a HC-Add-On make sure that the temperature sensor is inserted into one of the corresponding holes and close the door of your CompreDrive.



Approach your CompreDrive and proceed with your measurements.

# 11 Troubleshooting

## 11.1 No Connection to confocal system, system not found

Make sure the IP is correctly set in the confocal controller. Please also refer to the controller manual.

The correct IP has to be set in the connection dialog of RHD.MECServer.

Please also check the following points:

- Connect the sensor directly via the supplied USB-Ethernet dongle
- Ensure in the network settings, that the USB-Ethernet dongle has a static IP of 169.254.168.10 assigned
- Ensure that TCP/IP communication via port 1024 is not blocked by a firewall
- After plugging in the cable, make sure the two LEDs on the controller connector light up (green and yellow)  
→ If not, it means that the connection is not established correctly. Check the cable for damage.

## 11.2 Confocal system connected, but no data is transferred to Com-preDriveControl

Make sure the sensor is within its working range, and the light spot is visible on the mirror.

## 11.3 Data is present at the start point of a measurement, but is missing from some point

Make sure to use the full accessible measuring range. If the sensor is out of range due to sample compression, no data is transferred anymore.

## 12 Settlement

Warranty will be granted for a period of 2 years starting at the date of delivery.

Explicitly left out from warranty are parts that are subject to premature wear and tear due to use or other natural wear and tear. These components are regarded as consumables. Damage on the optical cable is excluded from warranty.

The costs for sending repaired or exchanged goods to the customer will be paid for by rhd instruments.

rhd instruments has to be notified of apparent defects and damages which occurred during production or delivery within 14 days after receiving the delivery. If a notification of apparent defects and damages does not occur within this period of time, the goods shall be deemed to have been accepted; as a result, the order will be assumed to be completed and approved.

Please note: Only workshops authorized by rhd instruments are allowed to perform repairs on the devices. If any mechanical or electronic components of the products are altered by customers themselves or by unauthorized workshops, a claim for warranty against rhd instruments is forfeited.

In case of a claim or sending back goods for repairs to be performed, please ask for the decontamination form beforehand. In general, rhd instruments must be contacted via e-mail or phone prior to any shipping of damaged goods.

## 13 Contact and Technical Support

For any questions with regard to our products, orders, or request for repairs please contact rhd instruments:

info@rhd-instruments.de

Phone: +49 6151 8707187

Fax: +49 6151 8707189

Web: <http://www.rhd-instruments.com>

rhd instruments GmbH & Co. KG

Otto-Hesse-Straße 19 / T3

64293 Darmstadt

Germany

Sitz der Gesellschaft: Darmstadt

Amtsgericht Darmstadt HRA 85824

WEEE-Reg.-Nr. DE 54715752

Haftende Gesellschafterin: rhd instruments Verwaltungs GmbH

(Sitz: Darmstadt, Amtsgericht Darmstadt HRB 96374)

Geschäftsführer: Dr. Benedikt Huber und Dr. Marcel Drüscher